



百特纳米粒度电位分析仪

体积分布报告页

V4.00

样品名称: M6563-30-2

样品来源: YP

检测人员:

检测单位: 河南省亚龙超硬材料有限公司

SOP名称: manual.sop

测试日期: 2024/3/26

测试时间: 13:22

备注: Z-均粒径为真实测量数值, 典型值为模拟数值, 仅供参考


分散介质: 水	样品材质: 聚苯乙烯乳球	检测角度: 90°	子测试次数: 14
温度: 25 °C	折射率: 1.59	平均光强: 363.60 kcps	衰减片编号: 4
粘度: 0.8936 mPa.s(cP)	吸收率: 0.01	分析模型: 通用	样品池种类: PS样品池
折射率: 1.33	激光波长: 671 nm	平衡时间: 120 s	检测点位置: 5 mm

测试结果	均径(d.nm)	面积(%)	峰标准差(nm)	CV (%)
Z-均粒径: 141.94 nm	峰1 104.52	100.00	37.51	35.89
PDI: 0.120	峰2			
截距: 0.79	峰3			

粒径(d.nm)	区间(%)	累积(%)
0.30	0.00	0.00
0.35	0.00	0.00
0.41	0.00	0.00
0.48	0.00	0.00
0.56	0.00	0.00
0.66	0.00	0.00
0.77	0.00	0.00
0.90	0.00	0.00
1.05	0.00	0.00
1.23	0.00	0.00
1.44	0.00	0.00
1.68	0.00	0.00
1.97	0.00	0.00
2.30	0.00	0.00
2.69	0.00	0.00
3.15	0.00	0.00
3.69	0.00	0.00
4.31	0.00	0.00
5.05	0.00	0.00
5.90	0.00	0.00
6.90	0.00	0.00
8.08	0.00	0.00
9.45	0.00	0.00
11.05	0.00	0.00

粒径(d.nm)	区间(%)	累积(%)
12.93	0.00	0.00
15.12	0.00	0.00
17.69	0.00	0.00
20.69	0.00	0.00
24.21	0.00	0.00
28.32	0.00	0.00
33.12	0.00	0.00
38.75	0.00	0.00
45.33	0.00	0.00
53.02	2.25	2.25
62.02	8.44	10.69
72.55	15.17	25.86
84.87	18.77	44.62
99.28	18.23	62.85
116.13	14.78	77.64
135.85	10.32	87.96
158.91	6.29	94.26
185.89	3.36	97.61
217.45	1.55	99.16
254.37	0.60	99.76
297.55	0.19	99.95
348.07	0.04	99.99
407.16	0.01	100.00
476.29	0.00	100.00

粒径(d.nm)	区间(%)	累积(%)
557.15	0.00	100.00
651.73	0.00	100.00
762.38	0.00	100.00
891.81	0.00	100.00
1043.22	0.00	100.00
1220.32	0.00	100.00
1427.50	0.00	100.00
1669.85	0.00	100.00
1953.34	0.00	100.00
2284.96	0.00	100.00
2672.88	0.00	100.00
3126.66	0.00	100.00
3657.48	0.00	100.00
4278.42	0.00	100.00
5004.77	0.00	100.00
5854.44	0.00	100.00
6848.35	0.00	100.00
8011.01	0.00	100.00
9371.05	0.00	100.00
10961.99	0.00	100.00
12823.02	0.00	100.00
15000.00	0.00	100.00

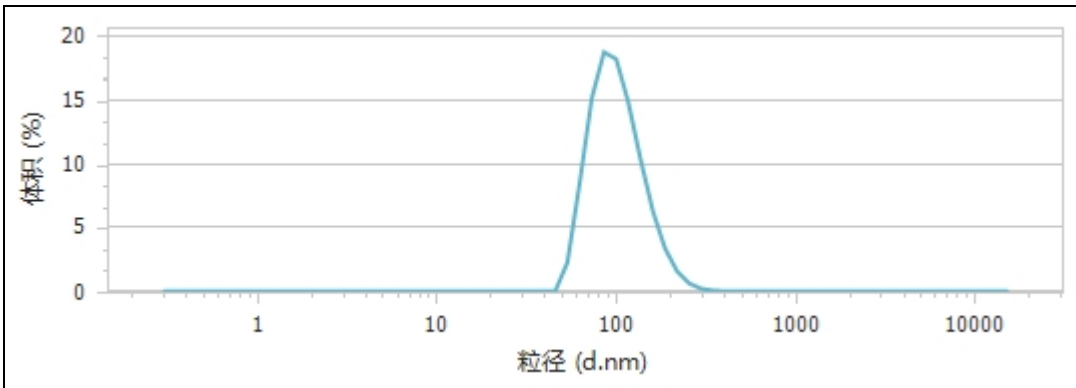

典型值	粒径(d.nm)
D10	61.29
D20	68.49
D30	75.27
D40	81.84
D50	89.12
D60	97.02
D70	107.43
D80	120.64
D90	143.32
D100	407.16