



百特纳米粒度电位分析仪

体积分布报告页

V4.00

样品名称: M6563-7-3

样品来源: YP

检测人员:

检测单位: 河南省亚龙超硬材料有限公司

SOP名称: manual.sop

测试日期: 2024/3/26

测试时间: 9:07

备注: Z-均粒径为真实测量数值, 典型值为模拟数值, 仅供参考


分散介质: 水	样品材质: 聚苯乙烯乳球	检测角度: 90 °	子测试次数: 19
温度: 25 °C	折射率: 1.59	平均光强: 286.07 kcps	衰减片编号: 4
粘度: 0.8936 mPa.s(cP)	吸收率: 0.01	分析模型: 通用	样品池种类: PS样品池
折射率: 1.33	激光波长: 671 nm	平衡时间: 120 s	检测点位置: 5 mm

测试结果	均径(d.nm)	面积(%)	峰标准差(nm)	CV (%)
Z-均粒径: 154.99 nm	峰1 122.02	100.00	35.70	29.25
PDI: 0.097	峰2			
截距: 0.85	峰3			

粒径(d.nm)	区间(%)	累积(%)
0.30	0.00	0.00
0.35	0.00	0.00
0.41	0.00	0.00
0.48	0.00	0.00
0.56	0.00	0.00
0.66	0.00	0.00
0.77	0.00	0.00
0.90	0.00	0.00
1.05	0.00	0.00
1.23	0.00	0.00
1.44	0.00	0.00
1.68	0.00	0.00
1.97	0.00	0.00
2.30	0.00	0.00
2.69	0.00	0.00
3.15	0.00	0.00
3.69	0.00	0.00
4.31	0.00	0.00
5.05	0.00	0.00
5.90	0.00	0.00
6.90	0.00	0.00
8.08	0.00	0.00
9.45	0.00	0.00
11.05	0.00	0.00

粒径(d.nm)	区间(%)	累积(%)
12.93	0.00	0.00
15.12	0.00	0.00
17.69	0.00	0.00
20.69	0.00	0.00
24.21	0.00	0.00
28.32	0.00	0.00
33.12	0.00	0.00
38.75	0.00	0.00
45.33	0.00	0.00
53.02	0.00	0.00
62.02	0.00	0.00
72.55	5.13	5.13
84.87	14.33	19.46
99.28	21.06	40.52
116.13	21.68	62.20
135.85	17.33	79.53
158.91	11.19	90.73
185.89	5.89	96.61
217.45	2.46	99.08
254.37	0.77	99.84
297.55	0.15	99.99
348.07	0.01	100.00
407.16	0.00	100.00
476.29	0.00	100.00

粒径(d.nm)	区间(%)	累积(%)
557.15	0.00	100.00
651.73	0.00	100.00
762.38	0.00	100.00
891.81	0.00	100.00
1043.22	0.00	100.00
1220.32	0.00	100.00
1427.50	0.00	100.00
1669.85	0.00	100.00
1953.34	0.00	100.00
2284.96	0.00	100.00
2672.88	0.00	100.00
3126.66	0.00	100.00
3657.48	0.00	100.00
4278.42	0.00	100.00
5004.77	0.00	100.00
5854.44	0.00	100.00
6848.35	0.00	100.00
8011.01	0.00	100.00
9371.05	0.00	100.00
10961.99	0.00	100.00
12823.02	0.00	100.00
15000.00	0.00	100.00

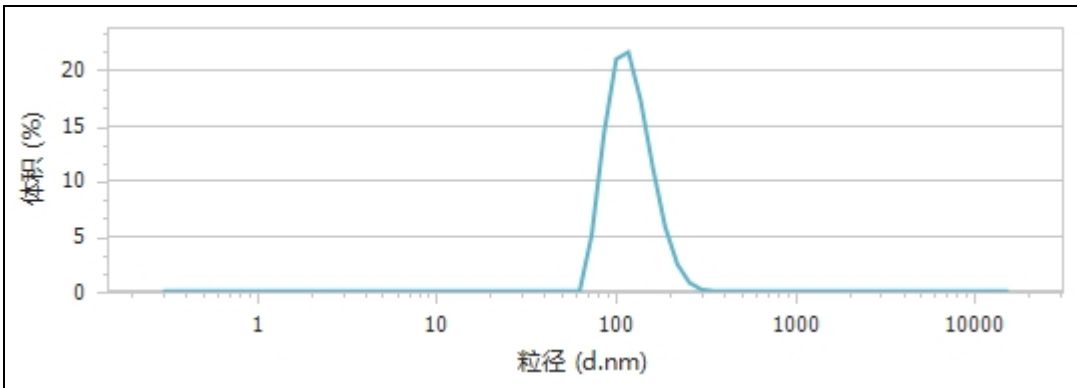

典型值	粒径(d.nm)
D10	76.74
D20	85.24
D30	92.08
D40	98.92
D50	106.65
D60	114.42
D70	125.00
D80	136.81
D90	157.41
D100	348.07